

Ambulatory Blood Pressure Monitor

TM-2430


Clinically
Validated



- 24 hour ambulatory blood pressure monitoring
- Light weight & compact size
- Oscillometric measuring method
- Battery Power (3 alkaline or Ni-Cd/Ni-MH "AA")
- Data output to computer or printer
- RS-232C Connectivity

Doctor Pro3 Blood Pressure Analysis Software

- Setting up all the parameters
- Various formats of graphs for further diagnoses
 - ◆ Trend
 - ◆ Correlation
 - ◆ Histograms
 - ◆ Circadian rhythm analysis
 - ◆ Hyperbolic index analysis

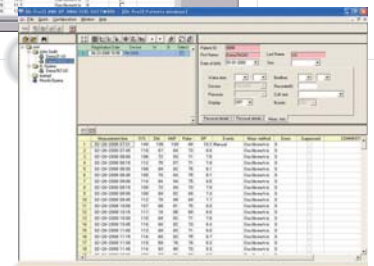

AND
A&D Medical
<http://www.aandd.jp>

Technology You Can Trust

Ambulatory Blood Pressure Monitor TM-2430

The TM-2430 is a compact & lightweight ambulatory blood pressure monitor that incorporates our latest technology gained from over 20 years of manufacturing ambulatory BPMs. This monitor was validated in accordance with British Hypertension Society protocol and accuracy was proven to be A/A. Weighing only 215 grams at the size of a cigarette pack, the TM-2430 can be carried by patients painlessly around-the-clock. Doctor Pro 3 Blood Pressure Analysis Software enables easy analysis of the data taken by the TM-2430, as well as parameter settings. (Optional)

- **24 hour ambulatory blood pressure monitoring**
- **Light weight & compact size**
- **Oscillometric measuring method**
- **Battery power (3 alkaline or Ni-Cd/Ni-MH "AA")**
- **Data output to computer or printer**
- **RS-232C connectivity**

- Automatically sets:
 - ◆ Proper cuff pressure
 - ◆ Exhaust speed
 - ◆ End of measurement
- New pulse correction algorithm for faster measurement & noise rejection
- Automatic or manual measurement
- Date/Time display: Until 2096 & auto leap year
- Multiple measurement intervals - adapts to patient's exact requirements.

Specifications

BP Measurement		
Measurement Method		Oscillometric
Measuring Range	Blood Pressure	40 ~ 280 mmHg (SYS - 60 ~ 280mmHg/DIA - 40 ~ 160mmHg)
	Pulse	30 ~ 200 bpm
Accuracy	Blood Pressure	±3 mmHg
	Pulse	±5 %
Pressurization	Method	Micro pump
	Maximum Pressure	300 mmHg
	Preset Pressure	Automatic
Depressurization	Method	Electric Controlled Exhaust Valve (ECEV)
	Safety System	A magnetic, solenoid valve opens when the pressure exceeds approximately 320 mmHg. Over 330 mmHg, an independent safety circuit exhausts air.
Interval Setting		5, 10, 15, 20, 30, 60, 120 min.
Clock		24 hour clock, calendar (until 2096)
Memory Data		Up to 350 sets of data
Power Supply	Battery	3x LR6("AA") type batteries - alkaline or Ni-cd type
	Battery Life	Approx. 200 measurements
Temperature & Humidity	Operating Temperatures	10°C ~ 40°C (50°F ~ 104°F), Less than 85%RH (non-condensing)
	Storage & Transportation	-20°C ~ 55°C (14°F ~ 131°F), Less than 95%RH (non-condensing)
External Dimensions		72(W) x 27(H) x 100(D) mm/2.8(W) x 1.0(H) x 3.9(D) inches
Weight		Approx. 215g (0.47lb) without cuff

Options and Accessories

- TM2430-02A Large cuff for left arm (28-36cm)
- TM2430-06A Adult cuff for left arm (20-31cm)
- TM2430-07A Small cuff for left arm (15-22cm)
- TM2430-09A Adult cuff for right arm (20-31cm)
- TM-2480 Printer terminal
- Doctor Pro 3 Blood pressure analysis software



TM-2480 Printer terminal



Doctor Pro 3

Blood Pressure Analysis Software (TM9501)

Doctor Pro 3 analysis software was designed to enhance diagnosing a patient's data taken by our ambulatory BP monitors (TM-2430) for further treatment. The following formats of graphs are available; Trend, Correlation, Histograms, Circadian Rhythm Analysis, and Hyperbolic Index (HBI or BP load) Analysis in addition to tabular data format. A file merge function is added to the latest version making closer analysis possible.

Also, the Doctor Pro 3 accepts UA-767PC data for long term diagnose.

TM-9501 Analysis software program

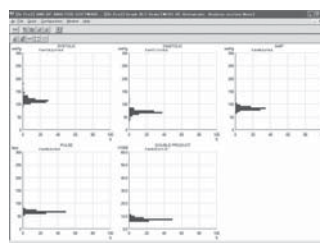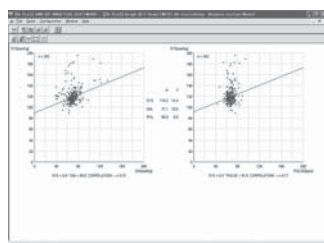

Doctor Pro 3 Hardware Requirements

CPU	Pentium III 1GHz equivalent or higher
RAM	256MB or more
Hard Disk	1GB required
Drive	CD-ROM drive
Interface	USB 2.0/1.1
OS	Windows™ 2000 Professional, Windows™ XP (Home & Professional) or Windows™ VISTA (all)

Windows™ is a registered trademark of Microsoft Corporation.



A&D Company, Limited
3-23-14 Higashi-Ikebukuro, Toshima-ku, Tokyo 170-0013 Japan
Telephone: [81](3)5391-6132 Fax: [81](3)5391-6148
1-243, Asahi, Kitamoto-shi, Saitama, 364-8585 Japan
Telephone: [81](0485)93-1111 Fax: [81](0485)93-1119
(EN46001 Certified Factory)
<http://www.aandd.jp>

A&D Medical/LifeSource
1756 Automation Parkway, San Jose, CA 95131 U.S.A.
Telephone: [1](408)263-5333 Fax: [1](408)263-0119

A&D INSTRUMENTS LTD.
Unit 24/26 Blacklands Way Abingdon Business Park,
Abingdon, Oxon OX14 1DY UNITED KINGDOM
Telephone: [44](1235) 550420 Fax: [44](1235) 550485
(Authorized Representative Established in the European Community)

A&D MERCURY PTY. LTD.
32 Dew Street, Thebarton, South Australia 5031 AUSTRALIA
Telephone: [61](8)8301-8100 Fax: [61](8)8352-7409

A&D RUS CO., LTD.
Verenskaya str.112 Kuntsevo Block, 121357,
Moscow, RUSSIA
Telephone: [7] (495) 937-33-44 Fax: [7] (495) 937-55-66

A&D KOREA Limited
Manhattan Bldg. 8F, 36-2 Yoido-dong, Youngdeungpo-gu, Seoul, KOREA
Telephone: [82](2) 780-4101 Fax: [82](2) 782-4280

A&D Instruments India Private Limited
509 Udyog Vihar Phase V Gurgaon-122 016, Haryana, INDIA
Telephone: [91](124) 471-5555 Fax: [91](124) 471-5599